



4
Linear
actuators



J X - 0 3
Gynecology Examination Table



Actuators With Min
IPX4 Protection
For Quiet Operation

Control Box
Built-in Battery



Calf Support



Remote Control
with Memory Function



Waste Basin

Soft, easy to clean
polyurethane
lying surface




Control units pictures are representative.



Foot Pressing

Leg Extension

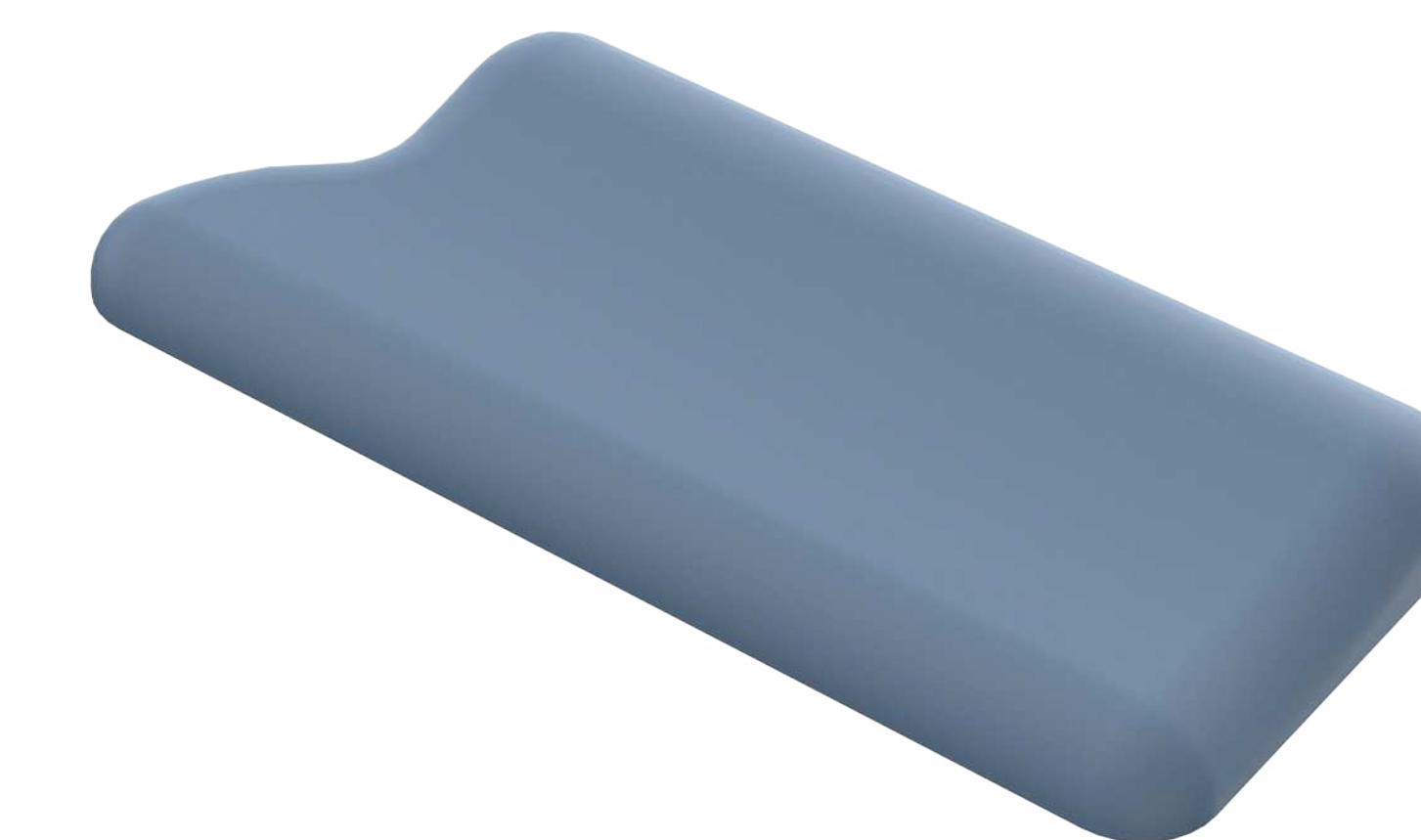

Arm Supporting



Height Adjustable Iv Pole



Optional Features



69-94 cm Height Movement



20° Trendelenburg / 3° Antitrendelenburg



52° Backrest Movement



Adjustable Calf support

Width (Over all)	Lenght (Over all)	Height Min.	Height Max.	Lying Surface	Sitting Surface	Backrest Angle	Trendelenburg Angle	Antitrendelenburg Angle	Safety Working Capacity	Operating Voltage	Number of Actuators
70 cm	131 cm	69 cm	94 cm	61x91 cm	58x36 cm	52°	12°	12°	230 kg	220 – 240 V AC	4

Soft, easy to clean
polyurethane
lying surface

ABS Plastic
Bottom Chassis
Cover

Nano ceramic coated
electrostatic powder
painted metal parts

Bottom chassis profiles:
30x50x2mm

Upper chassis profiles:
20x40x2mm

Anticorrosive
chromium-nickel
covered accessories.

Preventing wobbles
with a special lifting
column system reinforced
with a linear slide and carriage.

Ability to fixing
on the floor



Material Features



B. Kayacık Osb. Mahallesi
T. Ziyaeddin Akbulut Cad. Fabrika Üretim Apt.
No:19/A1 Selçuklu / Konya / Turkey
Tel : +90 332 239 10 80
Fax : +90 332 239 10 82 (pbx)
E-Mail : bilgi@medikal2000.com.tr
www.medikal2000.com.tr